

ageLOC™
boost
USER'S MANUAL

TABLE OF CONTENTS

CONTENTS	3
SAFETY INSTRUCTIONS, WARNINGS, AND PRECAUTIONS	4
BEFORE YOU BEGIN	6
CHARGING YOUR AGELOC BOOST DEVICE	6
USAGE INSTRUCTIONS	8
CLEANING AND CARING FOR YOUR AGELOC BOOST DEVICE	12
TROUBLESHOOTING	13
DEVICE DISPOSAL	13
MAINTENANCE	14
REPLACEMENT AND WARRANTY INFORMATION	15
PATENTS	15
DEFINITION OF SYMBOLS	15

CONTENTS

- 1 ageLOC Boost Device
- 1 Quick Start Guide
- 1 Magnetic Charger
- 1 Protective Conductive Surface Cover
- ageLOC Boost Activating Serum designed for use with ageLOC Boost device is sold separately.



SAFETY INSTRUCTIONS, WARNINGS, AND PRECAUTIONS

ageLOC Boost and Magnetic Charger are tested at pressures matching 4000 m in altitude.

This appliance contains batteries that are non-replaceable.

 **WARNING:** For the purposes of recharging the battery, only use the Magnetic Charger provided with this appliance.

-  The ageLOC Boost device is only for use on healthy skin. Do not use if you have open sores or wounds, are prone to rosacea, have allergies to metal, or if you have problematic or overly sensitive skin. Check with your physician before using if you are pregnant, have a pacemaker or similar device, are epileptic, wear metal braces on your teeth, have a metal implant, or in case of illness.
-  ageLOC Boost is designed for use by adults. Persons with reduced physical, sensory, or mental capabilities or lack of experience and knowledge and children above 16 years can use the device if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Keep out of reach of children. Children should be supervised to ensure that they do not play with the appliance. Cleaning and user maintenance should not be performed by children.
-  Do not apply directly to lips, eyes, or eyelids. The device and activating serum are for external use only. If you experience any swelling, prolonged redness or abnormal irritation, discontinue use immediately and seek medical advice if necessary. New ageLOC Boost users and those with sensitive or allergy-prone skin should always test on a small area of skin before using the product as directed.
-  To reduce the risk of shock, burn, fire, or injury, do not fully immerse your ageLOC Boost device in water. Do not place or store where it may fall or be pulled into a bathtub, shower, sink, or toilet.

-  Periodically inspect for damage; never operate the device if it appears to be damaged or broken to avoid injuries.
-  ageLOC Boost contains a lithium-polymer battery that is not replaceable. To prevent the possibility of explosion, never expose the device to heat. Do not store near a heat source, such as a radiator, fire, or heat vent. Check with your carrier before shipping or flying with this device.
-  Rechargeable batteries, like those of ageLOC Boost, need to carry a minimal charge in order to function properly. For optimal battery life, we recommend charging your device whenever low battery is indicated and to avoid storing the device unused for months at a time.
-  Note, ageLOC Boost is charged using inductive charging. Always use the provided charger when charging. Do not attempt to charge using any other inductive charger. If damaged, contact Nu Skin Support Services.
-  Do not charge this device in the bathroom or near water.
-  Do not use harsh chemicals or abrasives on your ageLOC Boost device.
-  Do not leave ageLOC Boost in an overly hot or cold environment for prolonged periods of time. Store in an indoor environment.
-  Always place your device on a heat-resistant, stable, flat surface, especially when charging.
-  Use ageLOC Boost only as directed and do not use on compromised skin or raised moles. Do not use ageLOC Boost for longer or more frequently than the application time indicated in the usage instructions. Regularly clean handgrip.

BEFORE YOU BEGIN


 Please read all safety and technical information before operating ageLOC Boost.

CHARGING YOUR AGELOC BOOST™ DEVICE

Fully charge ageLOC Boost before using for the first time.

Before charging, make sure both the device and the magnetic charger are completely dry.

Place the device on the magnetic charger (in the correct position, see diagram). Make sure the magnets snap into position. Plug the USB cable into a power source.

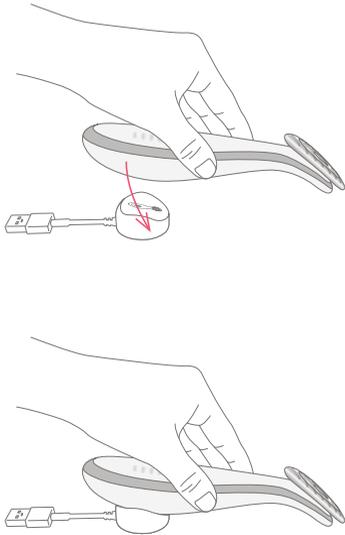

Charge ageLOC Boost completely.


When charging begins, ageLOC Boost will vibrate once and the LED light bar will slowly cycle up—this indicates the device is charging and displays the current level of battery charge.



After two minutes, the light bars will glow on and off to indicate the current battery level. Once the device is fully charged, the lights will stop glowing on and off and the LED light bars will remain dimly lit.



ageLOC Boost is charged using inductive charging. Always use the provided charger to charge the device. Do not attempt to charge using any other inductive charger. If damaged, contact Nu Skin Support Services.

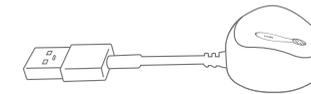


USAGE INSTRUCTIONS

Use ageLOC Boost™ as part of your morning skincare routine, after cleansing and toning.


PLEASE NOTE: ageLOC Boost has been specifically designed to be used with ageLOC Boost Activating Serum. For your safety and comfort, we do not recommend using the device with any product other than Activating Serum.

APPLICATION TIME


An application consists of two (2) one-minute, half-face application. After one minute, ageLOC Boost will pause briefly to allow you to apply additional ageLOC Boost Activating Serum to the conductive surface before you begin taking care of the second half of your face.

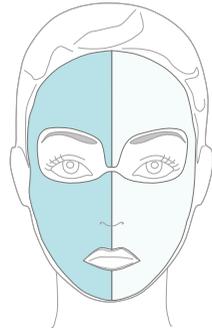


STEP 1 Turn on the device by swiping your finger up over the full length of the touch surface. Lights will turn on.

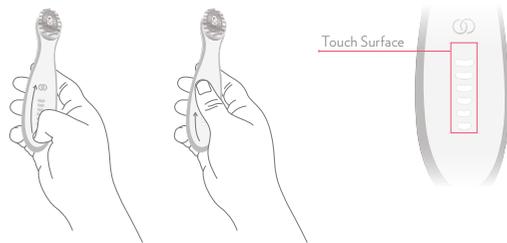


STEP 2 Dispense three pumps of ageLOC Boost Activating Serum onto the device's conductive surface. Using the device, spread Activating Serum evenly over one half of your face.

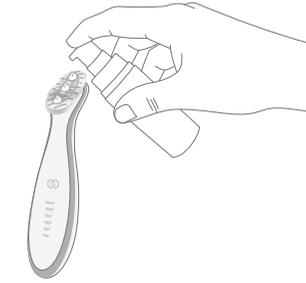

STEP 3 Begin application by swiping your finger up over the full length of the touch surface. The device will gently vibrate indicating application is in progress. Place conductive surface on face and begin the application on the skin by using slow movements in an upward and outward motion for one minute.

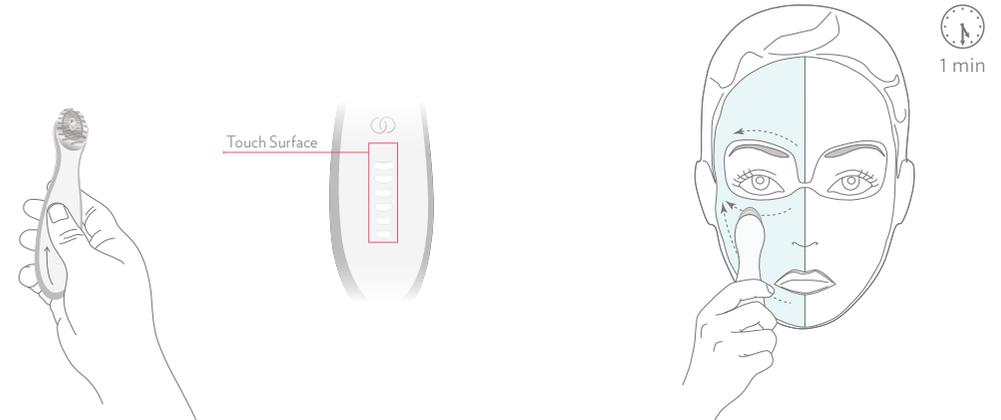


STEP 4 After the first one-minute application is complete, the device will pause automatically. Dispense an additional three pumps onto the conductive surface of the device, and spread evenly over the other half of the face. Swipe up the full length of the touch surface. The device will resume gently vibrating. Continue using slow, upward and outward motions for the remaining minute.

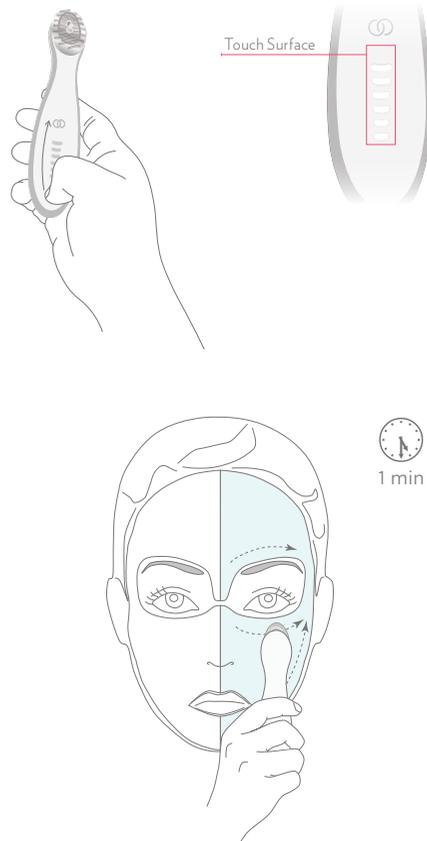

STEP 5 The device will automatically turn off after completing a full application cycle. Once application is complete, massage excess into skin, then follow with your daily skincare routine.

TO PAUSE APPLICATION At any time, simply swipe down over the full length of the touch surface. Swipe back up full length of touch surface to resume.

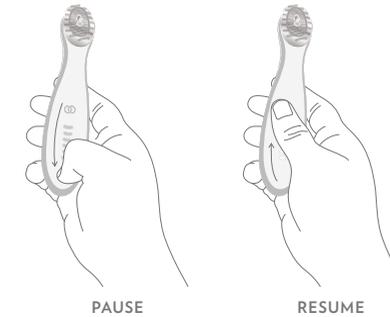


TO STOP APPLICATION MANUALLY During usage, you can stop the device at any time by swiping down the full length of the touch surface two times.



PLEASE NOTE: ageLOC Boost was designed to slowly glide across the skin's surface during use and should not be held in one spot. If you experience any prolonged redness or abnormal irritation, discontinue use immediately and seek medical advice if necessary.

CLEANING AND CARING FOR YOUR AGELOC BOOST DEVICE

ageLOC Boost should be cleaned according to these instructions after each use to avoid the accumulation of residue. After every application, rinse ageLOC Boost conductive surface under running water to ensure debris and product is removed from crevices. Thoroughly wipe the device dry with a soft washcloth. Store in a dry place.

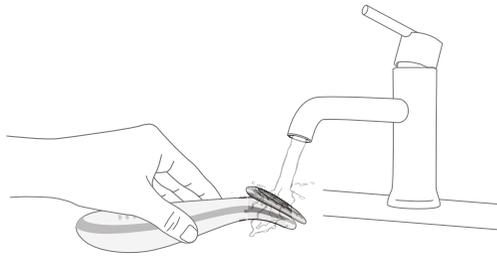

TROUBLESHOOTING

LOSS OF CONDUCTIVITY DURING USAGE: If the microcurrent conductivity is temporarily lost, the device's LED light bar will flash in continuous sets of five short flashes coupled with five sharp vibrations. If this happens, apply more ageLOC Boost Activating Serum to the device's conductive surface. Swipe up along the touch surface to resume application.

UNRESPONSIVE DEVICE: Swipe up slowly from the bottom of the device along the full length of the touch surface. If the device still does not respond, charge it fully and try again. If the device still fails to respond, contact Nu Skin Support Services.

LOW BATTERY: The bottom light on the device will turn amber. Charge the device.

EXTREME LOW BATTERY: The amber light at the bottom of ageLOC Boost will flash repeatedly when the user attempts to start an application. This indicates there is not enough battery remaining to perform an application, and the device will turn off. Charge the device.

DEVICE DISPOSAL

You must dispose of ageLOC Boost properly according to local laws and regulations. ageLOC Boost contains electronic components and a lithium-polymer battery, which requires disposal separate from household waste. When ageLOC Boost reaches its end of life, contact your local waste management department to learn about proper disposal and recycling options.



MAINTENANCE

Do not attempt to repair the device yourself. This will void any warranty. There are no user-serviceable parts.

Do not ship damaged devices by air.

TECHNICAL INFORMATION

ageLOC Boost

Maximum Output:

5 V, 250 μ A, 0.414 Hz

 IP25

 Charge only with ageLOC Boost Magnetic Charger.

If damaged contact Nu Skin Support Services.

ageLOC Boost Magnetic Charger

Input: 5 V  500 mA

 IP22

Electromagnetic Compatibility: Complies with the requirements of the EMC Directive 2014/30/EU

Safety: Complies with the requirements of the Low Voltage Directive 2014/35/EU

Weight of the device and battery pack: 69.5 g

MAINTENANCE, REPLACEMENT AND WARRANTY INFORMATION

LIMITED TWO-YEAR WARRANTY: Nu Skin® guarantees your device to be free from defects in materials and workmanship for a period of two years from the original date of delivery to you. This warranty does not cover damage to the product as a result of misuse or accident. If the product becomes defective within the two-year warranty period, please call your local Nu Skin® Support to arrange for a repair or replacement. You may be required to present a formal sales receipt showing the delivery date. Nu Skin® reserves the right to inspect the device. This warranty does not affect your statutory rights in law specific to the country of purchase, such rights remain protected. This warranty is valid in Austria, Belgium, Czech Republic, Denmark, Finland, France, Germany, Hungary, Iceland, Ireland, Israel, Italy, Luxembourg, Norway, Poland, Portugal, Romania, Russia, Slovakia, South Africa, Spain, Sweden, Switzerland, The Netherlands, Ukraine and United Kingdom.

PATENTS

Numerous U.S. and international patents pending.

DEFINITION OF SYMBOLS

Symbol	Definition
	Warning
	Consult User's Manual
	Dispose of Batteries Properly for Your Area
	Class III appliance supplied from a separated extra-low voltage (SELV) power source
	Charge only with specified charger
	Direct Current



CE EAC



NSE Products Europe BV, Da Vincilaan, 9 – 1930, Zaventem, Belgium

Imported by/Импортер/Импортер:

RU ООО "Нью Скин Энтерпрайзиз РС", 119049, Москва, ул. Донская 29/9, стр. 1, тел: 8-800-700-1984.

UA ТОВ "Ню Скін Ентерпрайзес Україна", 03680, Україна, Київ, вул. М. Грінченка, буд. 4, тел.: 0-800-30-19-84.

ZA Nu Skin Enterprises South Africa (Pty) Ltd. Techno Link 63, Regency drive, Route 21 Corporate Park, Nellmapius Drive, Irene, 0157, South Africa.

Made in China

www.nuskin.com

© 2021 NSE PRODUCTS, INC.

001931 Q1 2021